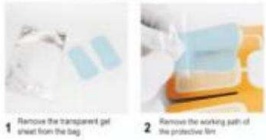
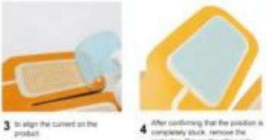
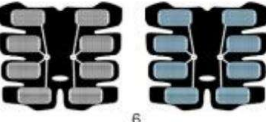
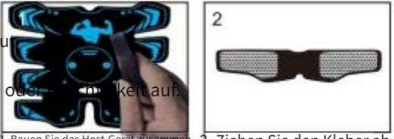
	HOW TO USE	HOW TO USE	HOW TO MANAGEMENT	
Directions for use  1 Remove the transparent gel sheet from the bag  2 Remove the working path of the protective film  3 to align the current on the product  4 After confirming that the position is completely stuck, remove the protective film on the other side   6	3-4. How to use an effective 1.The low frequency is suitable to 15 minuts one time and 2-3 times a day, if used it a long time or feel harder with impetus, then you will feel tired at moment. 2.It's able to avoid 30minutes after dinner, and favour to do after shower or warm washing. 3-5. Taking care after using 1.please turn off power after using. 2.please separate adhesive gel pads from machine. 3.Please attach adhesive gel pad to pad holder. (To prevent adhesive pad for drying, which must paste it pad holder) 4.The dirty point and etc on case should be cleaned by wet cloth by neutral cleaning material and wiped it by drying cloth again. 5.Don't keep in sunshine, hot temperature and humidity,wet and dust place, closer fire, or shocking place etc. 6.To keep it by children who can not touch. 7.Don't touch drying battery terminal with wet hand. 7	4-1. Main machine  1.The dirty point and etc on the case is cleaned by wet cloth by neutral cleaning material and wipe it by drying cloth again. 2.Do not keep it in high or low temperatures or humidity when storing the unit. 4-2. Adhesive gel pad 1.You will fell un-comfortable if the pasting side of adhesive pad gets dirty or pasting power is getting weakness more and more. Then, you can use wet towel and clean dirty side for regaining pasting again at moment, but, if it gets moisture too much, that will re-act to pasting. 2.The adhesive gel pads(electrodes) should be replaced with the same type and model. 4-3. Charging method 1. Plug the top of the main unit into the computer or mobile phone charger 2. The main unit light is on, indicating that it is charging, and the light is off,indicating that it is full. 8	4-4. Install instruction  1. Assemble the host device on abdominal pad  2. Peel off the adhesive pads.  3. Assemble the host device on arm pads  4. Attach the pads on your arms.  5.Stick the abdominal muscles onto the body 9	START FAT BURNING 

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